microlife

Sonic Electric Toothbrush

Clean S5

EN USER MANUAL

Please carefully read and consult the user manual before use and keep it for future reference.

Thank you for choosing Sonic Electric Toothbrush Microlife Clean S5 as your oral health care assistant.

Sonic Electric Toothbrush Microlife Clean S5 is an advanced electric toothbrush designed for effective and gentle cleaning of teeth, interdental spaces, gum massage polish of the teeth.

The sonic vibration motion technology allows you to reach the surface of the tooth even in the most hard-to-reach places.

The toothbrush has 5 modes:

- **Clean mode** – high speed oscillating toothbrush bristles can easily remove tartar and destroy hidden bacteria on the tooth surface and everywhere.
- **Sensitive mode** – gentle frequency targeting of sensitive teeth and gums.
- **Massage mode** – very gentle frequency of stimulation and massage of teeth and gums.
- **White mode** – switching between cleaning and polishing modes continuously helps to whiten and maintain the surface of teeth.
- **Polish mode** – ultra high frequency power swing to polish and brighten teeth.

Each mode works for 2 minutes. It is recommended to initially select one of the operating modes and remain in it until the end of the procedure

- Sonic Electric Toothbrush Microlife Clean S5 key functions:
- Wireless charging.
 - IPX7 water ingress protection.
 - Five operating modes for effective and gentle cleaning of teeth.

SAFETY INSTRUCTIONS

Please carefully read and consult the user manual before use and keep it for future reference.

WARNING!

- This device is equipped with built-in rechargeable battery. Do not throw it into fire, do not heat it up, and do not use or leave the device in places with high air temperature.
- Keep away from water or other liquids when charging.
- Keep the product away from strong electromagnetic field when charging.

- If any malfunctions are detected, stop using the device and contact a specialized service center. Do not attempt to repair the device yourself. All repairs must be carried out by an authorized service center. Unprofessional repairs may result in damage to the device, injury, and property damage.
- Do not disassemble the device.
- Only authorized service center personnel can replace the battery, unauthorized personnel are not allowed to disassemble the product.

- The manufacturer shall be not liable for any damage or loss caused by non-observance of safety instructions, improper or unintended use of the device.
- Do not use water higher than 40° C or acid cleaning solution to clean the brush handle, charging cable, power cord and adapter, the adapter is strictly prohibited in water.
- Before cleaning the charging device, be sure to pull out the adapter from the socket, otherwise it may cause electric shock.


CAUTION!

- The toothbrush is not a toy and is not intended for play. Sonic Electric Toothbrush Microlife Clean S5 contains small parts which can be swallowed by children. For safety reasons, this device may not be used by children, persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children should not use this product without adult supervision. Keep the product and cable out of the reach of children.
- Please check the product for obvious damage before each use. If the toothbrush is damaged or seriously dropped, please stop using it.
- Do not charge in the bathroom or humid environment, and do not charge this product for a long time to avoid unnecessary dangers.
- Basic precautions must always be observed when using any electrical device.

PRECAUTION!

- This device is intended only for household use. Do not use the device for purposes other than intended. Use the device only for the purposes described in this User Manual. The use of the device for purposes other than those specified in this manual is considered as a violation of the operation rules.
- During the first days of using the toothbrush, you may feel irritation or tingling in the tooth area, and you may experience some bleeding or discomfort in the gum area due to the high-speed vibration of the bristles. In general, bleeding should stop after a few days. Should it persist after two weeks, please consult your dentist or dental hygienist.
- Consult your dentist before using a toothbrush if any of the following concerns:
 - sensitive teeth or gums;
 - recent periodontal or dental surgery;
 - undergoing treatment for oral problems.
- Stop using the product or consult your dentist if you experience any discomfort or pain while using the sonic toothbrush.
- If your toothpaste contains peroxide, sodium bicarbonate, or other bicarbonate (a common ingredient in whitening toothpaste), be sure to clean the brush head with water after each use to prevent any plastic breakage.
- The toothbrush is a personal hygiene product and must be used individually.
- Do not press the brush too hard on your teeth and gums.
- The product is intended for indoor use only. Keep it out of direct sunlight and store it at room temperature or in a dry and ventilated place.
- Never place or insert any object into the charging port of the product.
- Do not connect the charger to the mains with wet hands.
- Only connect the charger to an electrical mains supply with the required parameters specified in this manual.
- Switch off the device after each use.
- The toothbrush head should be replaced every 3 months or more frequently when the bristles become worn out.
- Use only brush heads recommended by the manufacturer.

DEVICE DESCRIPTION

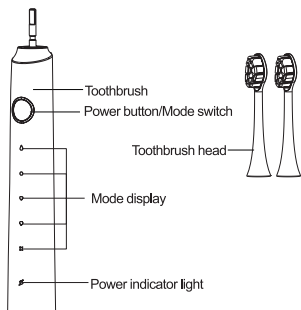COMPLETE SET

Sonic Electric Toothbrush Microlife Clean S5 (handle) – 1 pc.
Replaceable brush head* – 2 pcs
(ultra-thin – 1 pc, medium – 1 pc).
Charging base – 1 pc.
User manual – 1 pc.
Storage – 1 pcs.

*Ultra-thin brush head consist of specially designed ultra-thin bristles less than 0.01 mm for superior cleaning of hard-to-reach areas, inc. between teeth, efficiently sweeping away food and plaque. Medium brush head consist of high-density bristles with tap rounding treatment brush wire. This ensures thorough cleaning without harming the gums.

The toothbrush head should be replaced every 3 months or more frequently when the bristles become worn out.

BASIC SPECIFICATION

Type	Sonic Electric Toothbrush
Model	Microlife Clean S5
Dimensions (L x W x H), mm	245 x 26 x 26
Weight, g	86
Maximum vibration frequency (movements per minute)	~42000
Rated Voltage / Power (Charging Base)	USB 5.0V --- 3.0W
Battery type / Power source	Battery IMR 14500, 500 mAh 3.7 V
Charging Time	Wireless charging approximately 9 hours
Waterproof rating	IPX7
Storage conditions*	From +5°C to +45°C with relative humidity up to 65%
Transport conditions*	From -10°C to +50°C with relative humidity up to 75%
Operation conditions	From +15°C to +35°C with relative humidity up to 65%
Stiffness bristles	Ultra-Soft Bristles, Medium Bristles

WARNING!


*After transportation or storage at low temperatures, it is necessary to keep the device at room temperature for at least 2 hours before switching on.

STARTING TO USE

Read the user manual carefully before using the device and keep it for future reference.

- This product has a one-button start system, including ON/OFF, mode switching.
- 2 minutes intelligent timer: after using any mode for 2 minutes, the toothbrush automatically turns OFF.
- 30 seconds zone change reminder: pause 0.3s every 30 seconds for each mode.
- To switch modes, press and hold the Power Button pressed for 5 seconds to switch modes. After 5 seconds from the moment the toothbrush starts working, press the Power Button if you want to turn off the device.


1. Align the front part of the toothbrush head with the front part of the toothbrush handle.
2. Install the toothbrush head on the metal shaft (there will be a small gap between the two components, to prevent the brush head from colliding with the main engine during work, so that the motor has more powerful power).

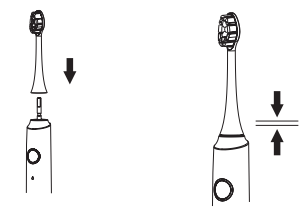3. Dampen the bristles and apply some toothpaste to the bristles.

4. Gently press the Power Button to turn ON the toothbrush. Place the toothbrush bristles against gum line at a slight angle (45 degrees), moving gently in circular motion from gum line, pressing gently to make the bristles reach the area between beneath the gumline and teeth.

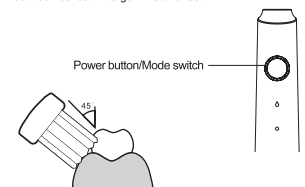This product has an intelligent memory function, which automatically enters the previous usage mode when activated next time.


The toothbrush has 5 modes:

- **Clean mode** – high speed oscillating toothbrush bristles can easily remove tartar and destroy hidden bacteria on the tooth surface and everywhere.
- **Sensitive mode** – gentle frequency targeting of sensitive teeth and gums.
- **Massage mode** – very gentle frequency of stimulation and massage of teeth and gums.
- **White mode** – switching between cleaning and polishing modes continuously helps to whiten and maintain the surface of teeth, making teeth as white as porcelain.
- **Polish mode** – ultra high frequency power swing to polish and brighten teeth.

Each mode works for 2 minutes. To switch modes, use the Power Button.

5. Brush your teeth with small back and forth motion so the bristles reach between the teeth, continue this motion throughout your brushing cycle. Move up and down, do not press hard or scrub.

6. Start brushing from area 1, then move to the next area so the bristles reach between the teeth, continue this motion according to dentists' recommendations. After every 30 seconds, the brush takes a short pause of 1–2 seconds, thus reminding you to change the cleaning area. After the 2nd minute, the brush turns off.

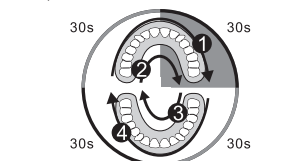7. Rinse the toothbrush head under a stream of water for a few seconds after each use. After cleaning the brush head, place the toothbrush in dry, cool and ventilated, dry place. If necessary, wipe all surfaces with a soft cloth. Do not use sharp objects for cleaning.

WARNING!

- Do not put the heads and the handle in the dishwasher.
- Do not wash the device in boiling water.

CHARGING DESCRIPTION

1. Please check the power supply voltage and the power socket are safe and normal and the voltage adaptation.
2. Please make sure that the USB cable is dry, the charging plug is free from water stains and damage and then plug power supply under safe conditions.

3. When using an electric toothbrush for the first time, please turn on the electric toothbrush first. In case of low or no power, please fully charge the electric toothbrush first.

BATTERY INDICATOR DESCRIPTION

Indicator status during discharge	Description
Indicator ON	Standby Modes
Red indicator ON	Low battery, please charge in time
Indicator not ON	Ultra low battery, products cannot work normally

Indicator status during charging	Description
Red indicator ON	On charging
Indicator turns from red to green	Finished charging

WARNING!

- Please pay attention to the status of the battery indicator during use.
- Please charge in time when the battery is low.

If you have not used the electric toothbrush for a long time, please charge the toothbrush when you use it again to ensure safe charging, which is conducive to improving the service life of the product.

If the product is not charged for a long time, the indicator light may not immediately light up during charging, which is a normal situation. Generally, it needs to be pre charged for 30 minutes before the indicator light can work properly.

TROUBLESHOOTING

Problem description	Possible cause	Solution
The toothbrush does not turn on	Low battery power	Charge the toothbrush
Weak vibration of the brush head	Low battery power	Charge the toothbrush

DISPOSAL

- This symbol on the product or its package means that this product does not fall under the category of domestic waste.
- Proper disposal of the device will prevent adverse environmental and human health effects.
- To protect the environment, the device must not be disposed of together with domestic (household) waste. Disposal shall be provided in accordance with local regulations.
- The device must be disposed of in accordance with the EU Directive 2012/19/EU WEEE (Waste Electrical and Electronic Equipment).

If you have any questions, please contact the local public utility responsible for waste disposal.

APPLIED STANDARD

Sonic Electric Toothbrush Microlife Clean S5 complies with:

- Directive 2014/30/EU EMC;
- Directive 2011/65/EU RoHS with its Amendments (EU) 2015/863.

LIST OF SYMBOLS

Symbol	Description
CE	CE marking, Conformité Européenne (European Conformity). This symbol means that the device fully complies with applicable European Union Acts and European technical conformity
WEEE	WEEE (Waste Electrical and Electronic Equipment Directive). The symbol on the product or its package means that this product does not fall under the category of domestic waste. To avoid possible harm to the environment and human health, separate such wastes from others and dispose of them in accordance with accepted standards.
Manufacturer	Manufacturer
Date of manufacture	Date of manufacture
Consult User Manual before use	Consult User Manual before use
Class II equipment. Protection against electric shock	Class II equipment. Protection against electric shock
Not suitable for children under 3 years. The device is not a toy and is not intended for games	Not suitable for children under 3 years. The device is not a toy and is not intended for games
Ingress protection degree provided by IEC 60529	Ingress protection degree provided by IEC 60529
IPX7	IPX7
Indicates the entity importing the medical device into the locale	Indicates the entity importing the medical device into the locale

WARRANTY INFORMATION

The warranty period of this product is 2 years from the date of sale, providing that all operation conditions specified in this User Manual are strictly observed. The warranty is valid only on presentation of the sales receipt and the warranty card completed by the authorized dealer. The warranty does not apply to consumables, constituent parts and components, such as battery, package, etc.

Warranty and free maintenance service are not provided in case of:

- violation of User Manual.
- damage caused by intentional or erroneous actions of the consumer due to improper or negligent handling;
- evidence of mechanical impact, dents, cracks, chips, etc. on the device housing, evidence of opening the device, disassembly, unauthorized repair, ingress of moisture, exposure to aggressive agents or any other unauthorized alterations of the device, and in other cases of violation of storage, cleaning, transportation and operation rules specified in this User Manual.
- ingress of oils, dust, insects, liquids and other foreign objects inside the device.

Follow the instructions precisely to ensure reliable and long-term operation of the device.

In case of abnormal operation of the device, please contact your retailer or Customer Service.

The manufacturer reserves the right to make alterations in the design of the device.

WARNING!

- Do not try to repair the device yourself – this will void the warranty.
- The device alterations without the manufacturer's permission will void the warranty.
- For repair and maintenance, please contact a specialized after-sales service.

Manufacturer: Shenzhen Samo Technology Co., Ltd., 2F Building No. 3, He Sheng Park, He Ping District, Fu Hai Street Bao'an, Shenzhen, 518000 P.R.C

Importer in the EU: Samocontrol Ltd, 2/6 Konstantin Velichkov St., Plovdiv 4000, Bulgaria

Microlife is a registered trademark of Microlife AG, Switzerland.

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Sonic Electric Toothbrush Microlife Clean S5

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The set completeness and operation of the device are tested in the presence of the buyer.

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ATTENTION! This warranty is valid only if documents are filled in properly.

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